




The STERIS 5085 Harmony Surgical Table is a mobile, electrohydraulically operated surgical table designed to support virtually all general and bariatric surgical procedures including cardiac and vascular, endoscopic, gynecology, urology, nephrectomy, neurology, ophthalmologic and orthopedic with the addition of STERIS table accessories.

Designed to quickly and easily adapt to virtually every surgical specialty, the 5085 Harmony Surgical Table offers a higher weight limit, exceptionally low height and a large imaging area.

The surgical table supports infection control measures with a hook and loop fastener-free table pad attachment and flush rail design that eliminates exposed screws, making it easier to clean.

Features

- Accommodates virtually all patients with generous weight and height ranges
- Provides unrestricted perineal access with removable leg section
- Optional Featherweight Leg Section complies with AORN Safe Handling principles
- Enhances lateral procedures with a powered, radiolucent kidney elevator
- Provides outstanding access for C-arm, permitting clear high-quality images crucial to MIS, cardiothoracic, orthopedic and neurosurgical procedures
- Facilitates quick and safe patient positioning with an easy-to-use, easy-to-read hand control
- Contours to virtually all patient postures with four section tabletop
- Compensates for uneven OR floors with self-leveling floor locks
- Prevents sectional and table base collisions with the Auto Limit Sensor
- Compatible with a wide array of previously purchased accessories

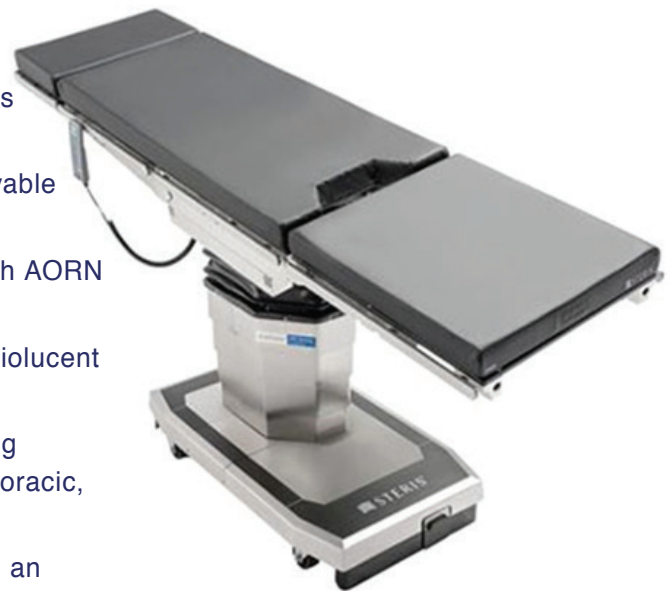

Specifications

Technical Details

Table Length: 81 in (2057 mm)
Table Width: 20 in (508 mm) (28° with extenders)
Table Height Range: 22.8 in to 42.5 in (581 to 1080 mm)
Table Slide Range: 9 in (227 mm) to head, 9 in to foot (18 inches total)
Patient Weight Capacity:
1,200 lbs. (544 kg) patient support, including raise/lower (centered on the column)
900 lbs (454 kg) full table articulation (centered on the column)
600 lbs (272 kg) full table articulation, including slide
Table shipping Weight: 820 lbs. (372 kg)
trendelenburg/Reverse: 30°/30°
Lateral Tilt(left/Right): 20 °/ 20°
Head Section: +90° to -90° (15° increments)
Back section: +80° to -40°
Leg Section: 0/-105° (removable)
Flex/Reflex: 140°/100°
manual Override: yes
Perineal Cut-out: yes
Kidney Elevator: 4in (110 mm)powered