

MEDIVATORS ADVANTAGE PLUS




ABOUT

The **Medivators Advantage Plus** makes reprocessing flexible endoscopes easier and more reliable. This fully automated system takes care of high-level disinfection while following strict infection control standards to help protect patients from cross-contamination. Its dual-basin design allows two endoscopes to be processed at the same time, improving turnaround and reducing downtime. The system also includes built-in software that tracks and records each cycle, simplifying compliance and reporting for your quality team.

FEATURES

- **Dual-Basin Efficiency:** Reprocesses two endoscopes at once to save time and increase productivity.
- **Automated HLD Cycle:** Provides a consistent, validated disinfection process that meets AAMI, AORN, and CDC requirements.
- **Integrated Data Management:** Records, stores, and exports cycle information for compliance and quality assurance.
- **User-Friendly Interface:** Color touchscreen makes operation simple, with clear cycle status updates.
- **Advanced Monitoring:** Tracks key parameters like chemical concentration, temperature, and pressure for reliable results.
- **LEAK TESTER 9000™ Compatibility:** Optional integrated leak testing helps extend the life of your endoscopes.
- **SQL Database:** Securely stores cycle data for easy access and audit readiness.

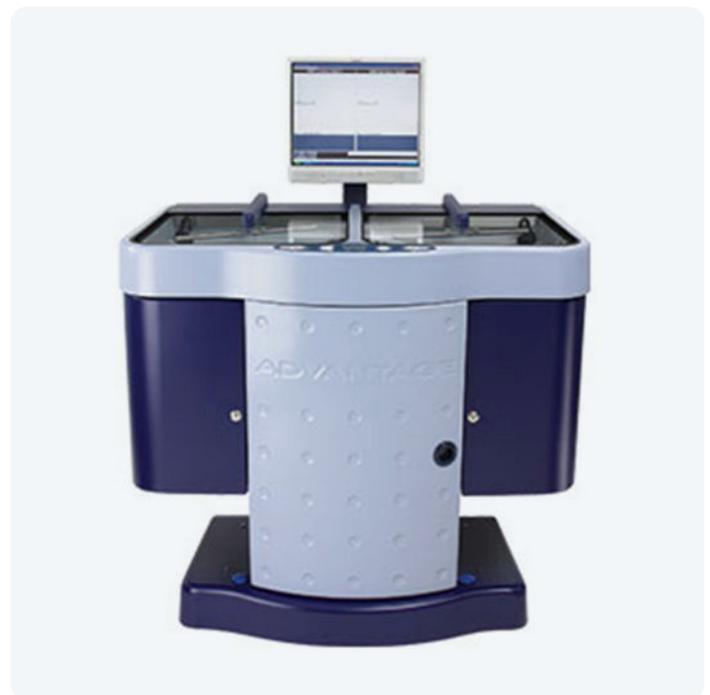

SPECIFICATIONS

DIMENSIONS

Height: 49.5 in (125.7 cm)

Width: 28.5 in (72.4 cm)

Depth: 29.5 in (74.9 cm)

Weight: 385 lbs (174.6 kg)

POWER

Voltage: 100-240V AC

Frequency: 50/60 Hz

Input Current: 10A (115VAC)

ENVIRONMENTAL CONDITIONS

Operating Temperature: 60°F to 90°F (15°C to 32°C)

Relative Humidity: 20% to 80% (non-condensing)

Water Pressure: 20-80 psi (138-552 kPa)

Water Temperature: 60°F to 85°F (16°C to 29°C)

PERFORMANCE

Water Quality: Requires bacteria-free, filtered water (0.2 micron filter)

Drain Connection: Requires a dedicated, air-gapped drain

Air Supply: Requires medical-grade, oil-free, dry compressed air