

# HILLROM LIFT-ASSIST BEACH CHAIR




## ABOUT

The Hillrom Lift-Assist Beach Chair is a surgical table accessory used to position patients in a modified Fowler posture for shoulder and upper body operations. It provides firm support for the head, torso, and shoulders, giving surgeons clear access to the operative site. Built with caregiver safety in mind, the chair reduces the effort needed for manual positioning. It fits most international operating table rail systems, allowing easy integration into orthopedic and general surgical environments.

## FEATURES

- **Lift-Assist:** A single release handle controls the motion for raising or lowering the patient, enabling smooth adjustments with minimal strain on staff.
- **Mounting System:** Specialized clamps fasten tightly to US, EU, UK, and JP table rails,
- **Shoulder Panels:** Lockable panels give lateral support and can be removed from the operative side when unobstructed access is required.
- **Head Support:** An extendable head positioner and tube assembly keep the head and neck secure as the patient's position changes.
- **Weight Capacity:** Supports up to 500 lbs (226 kg), making it suitable for a wide range of patients.
- **Restraint Compatibility:** Works with traction straps and arm supports to stabilize both the patient and the non-operative limb during surgery.

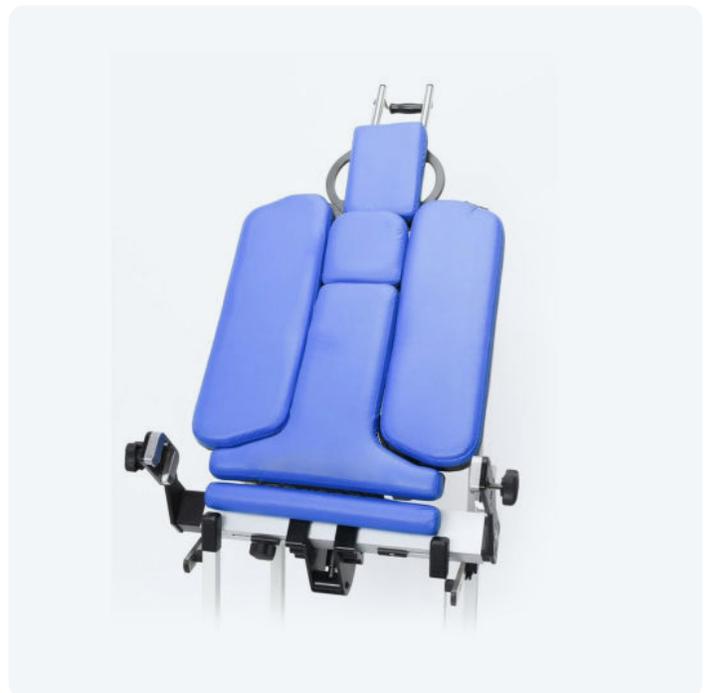

# SPECIFICATIONS



## DIMENSIONS

- **Width:** 26 in (66 cm)
- **Height:** 13.75 in (35 cm)
- **Length:** 38 in (96.5 cm)
- **Overall Weight:** 35 lbs (15.9 kg)
- **Safe Working Load:** 500 lbs (226 kg)

## ENVIRONMENTAL CONDITIONS

**Storage Temperature:** -29° C to +60° C

**Storage Relative Humidity:** 15% to 85%

**Operating Environment:** Controlled Operating Room conditions.

## COMPATIBILITY

**Compatible Rail Styles:** US (3/8in x 1 1/8in), EU (10mm x 25mm), UK (6.35mm x 31.75mm), DEN (6.35mm x 38.1mm), JP (5.9mm x 32mm), SWISS.

**Key Accessories:** Lift-Assist Beach Chair Clamps (A-91503), Intubation Plate and Pad (A-91502), LPS Arm Support (A-90002).

## MATERIALS

**Primary Materials:** Stainless Steel, Aluminum, High-Density Foam, Fiberboard, Delrin, UHMW Polyethylene, Nylon, Polyurethane.